


Chamlide PM

Perfusion type multi-well closed magnetic chamber



- Uses a 24 x 60mm coverslip as a bottom for high magnification microscopy.
- Closed type but perfusion is also possible.
- Can individual 5experiments on a coverslip by the media perfusion.
- Patented magnetic system for easy assembly.
- Control the temperature of the chamber using a heating plate.
- Use with the fluidic inline heater to control medium temperature.
- The height of the channel can be adjustable.

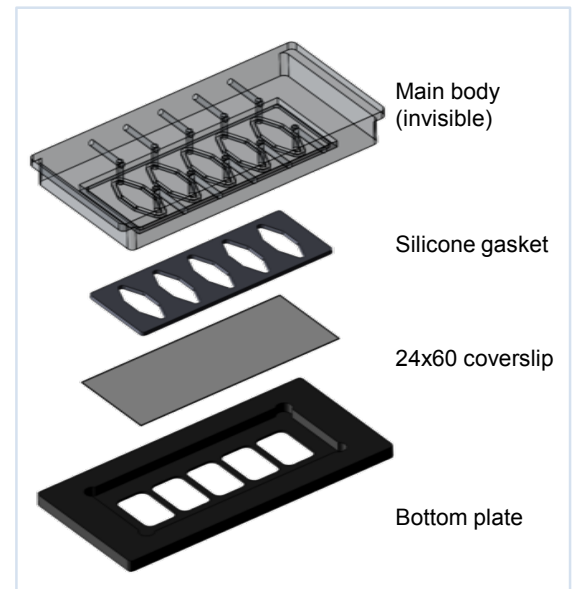

Specifications

Magnetic closed chamber internal volume

Depends on the channel size

Material

Main body	Acryl
Bottom	Black anodized aluminum alloy
Gasket	Silicone

Model No.

PM-S-10

Product

Chamlide PM