
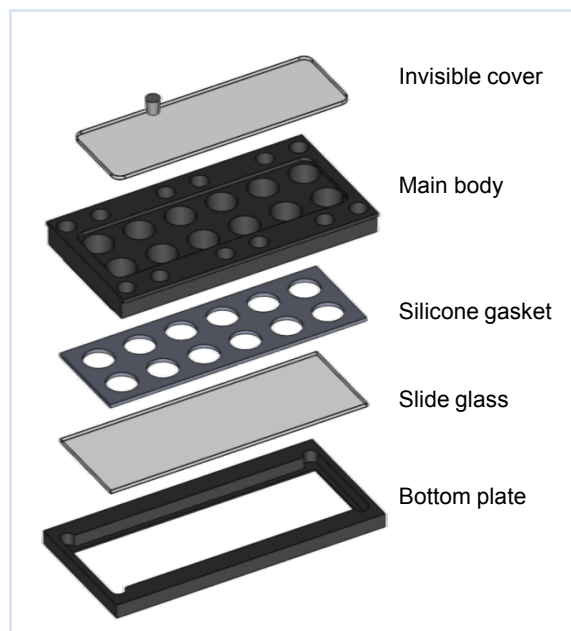


# Well-Slide

## Slide glass bottom multi-well magnetic chamber

- Slide glass bottom multi-well chamber
- Best for immunohistochemistry staining.
- Replaces an expensive consumable chamber slide.
- Semi-permanent (re-useable).
- Simply remove and replace the bottom slide glass for next use.
- Autoclavable.
- Uses silicone gaskets to prevent leaking.
- Metal bottom for stable temperature control during live cell imaging.
- Non-toxic, non-reactive PEEK material.
- Patented magnetic system for easy assembly.
- Control the temperature and the CO<sub>2</sub> level using Chamlide TC for live cell imaging.

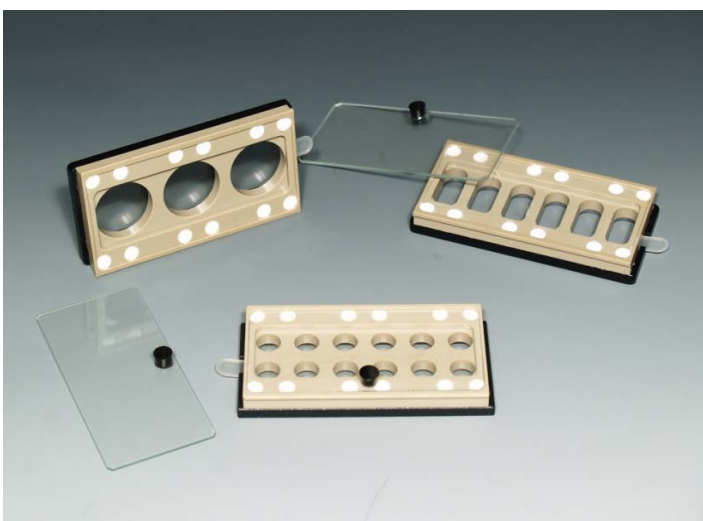

### Specifications

Well-slide internal volume	12 well type	400 $\mu$ l / well
	6 well type	1 ml / well
	3 well type	2 ml / well
Material	Main body	PEEK
	Bottom	Black anodized aluminum alloy
	Gasket	Silicone

Model No.	Product
WS-S-03	3 Well-slide
WS-S-06	6 Well-slide
WS-S-12	12 Well-slide